

Musician ER Earmolds



Oto Seal



- Beige (default)
- Available in all colors
- Clear if requested
- Filters available 9, 15, 25
- Countersunk if requested

Silicone



- Opaque in color (default)
- Available in all colors
- Countersunk only


Musician ER Earmolds

High-Fidelity Hearing Protection

Music • Clubs • Concerts

OTO Hearing Products is an authorized manufacturer by Etymotic Research to produce ER Musician ear plugs. Musician plugs are required to use an acoustic mass meter to verify that these plugs have the correct volume of air in the sound bore, which produces flat attenuation. Here at Oto, we manufacture these in two different materials, one being OTO Seal which is a rigid material that will become semi-flexible from body temperature when worn, and the other from a soft high grade silicone.

Musicians practice and perform in a variety of different settings and they are exposed to high levels of sound, sometimes for long periods of time. They require different amounts of protection, depending on the sound levels they encounter during rehearsals and performances. As designated by the number classification, ER 25's and ER 9's provide either more or less attenuations than ER 15's. Some musicians use one type of attenuator in one ear and one type in the other ear, depending on the source and location of the sound.

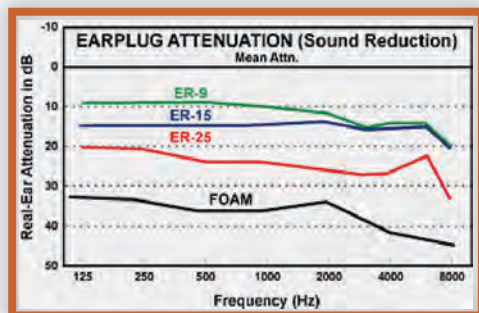

- Quality is clear and natural
- Sound not muffled
- Fatigue from noise is reduced
- Filters are interchangeable
- Timbre and richness of music are preserved

www.otohearing.com

Oto Hearing Products
1-800-661-7723 Toll Free
604-590-1211 Local